



independant® Care
Healthcare Design Technology

Why choosing **independant®** Care?

Section Configuration

Green Lighting

Decrease your electricity consumption thanks to our innovative LED lighting strips which are available for all direct, indirect and even night lighting. Lower maintenance, greater performance! Up to 40% energy saving; 50,000 hours average lifetime.

Powder coating, powered hygiene

Epoxy powder coating guarantees a clean smooth hygienic surface for all our units. Our dedicated antibacterial paint can be added to the powder coating as an option to ensure a high level of safety against bacteria.

Easy to install easier to maintain

We create our equipment based on two simple key principles: ease of installation and long-lasting top quality components ensuring minimum maintenance and maximum reliability. independant® Care is delivered fully assembled to reduce costs of installation.


Modular versatility

independant® Care is arranged to suit perfectly any medical requirement. The modularity of the unit allows any specification to be met with one single product which can be adapted to any level of care needed in the hospital.




Version ES

This is the starting version of independant® Care available in horizontal or vertical orientation. Ideal for all standard wards, it offers a simple design featuring a single profile to house both electrical and gas components whilst maintaining segregation between the two. The ES range can be offered with integral lighting, additional DIN rail and as a horizontal or vertical orientated system.

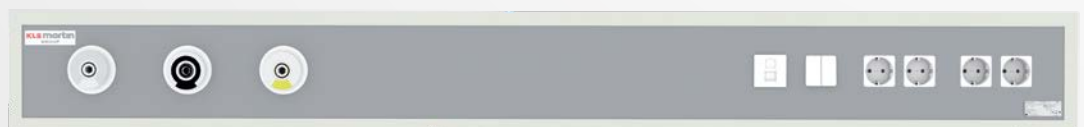


Version MD

Suitable for both standard and critical care this version features two separate profiles housing all services. Electrical and medical gas components run in completely separate compartments and are independently fed to ensure segregation is maintained. Similar to the ES model, it can be fitted with lighting modules, DIN rails and depending to the configuration in horizontal or vertical orientation.

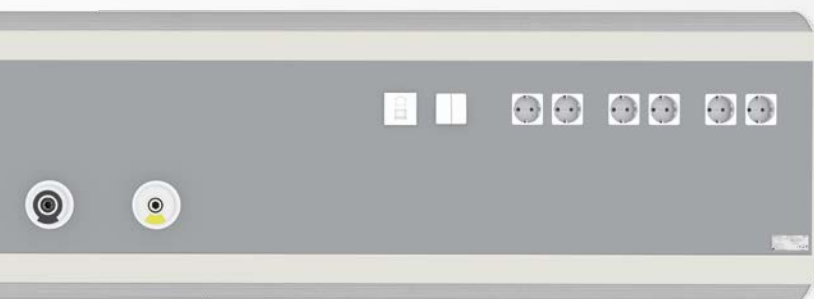


Version LG

Ideal for all critical areas this configuration includes three profiles to guarantee that even with large number of services segregation can always be maintained. Although generally supplied without lights both direct and indirect lighting can be included as option together with DIN rails if needed. According to your needs the system is in a horizontal or vertical orientation available.

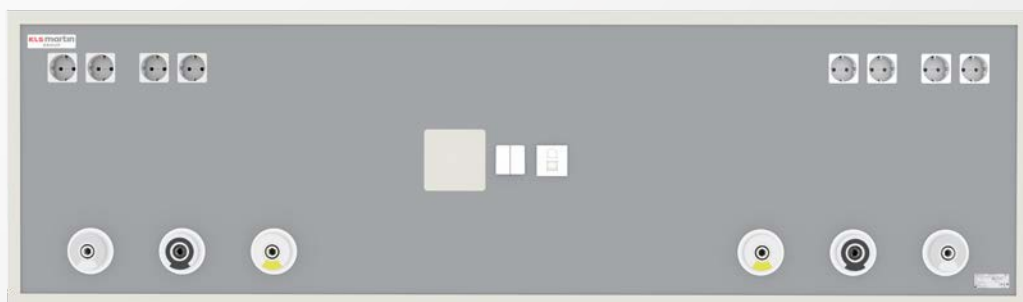

Horizontal wall mounted

The profile overview explains the various configurations available for each independant® Care Horizontal model. In all profiles medical gas components are kept segregated from the electrical outlets.

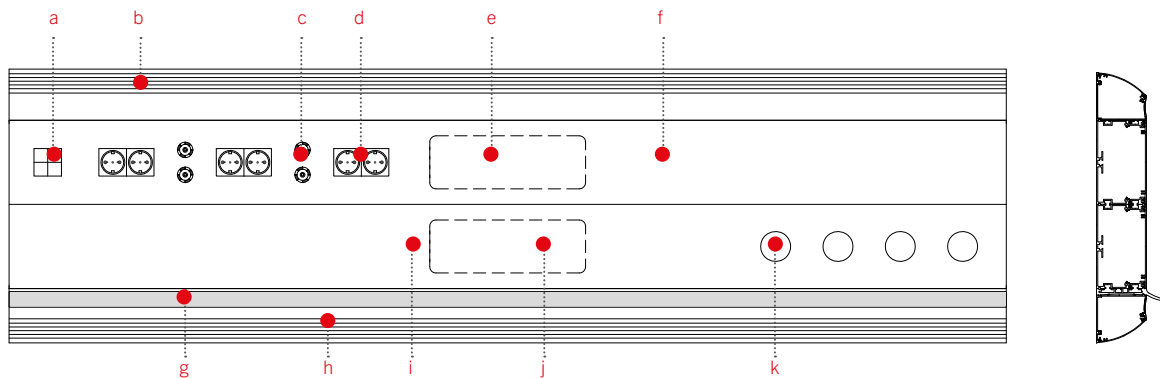


This image is for general demonstration only. The number of services varies according to the specifications given by each customer.

- | | |
|---|----------------------------------|
| a Switches | g DIN rail |
| b Indirect light | h Direct light |
| c Equipotential ground socket | i Compartment for gas components |
| d Electrical socket | j Gas entry fed |
| e Electrical entry fed | k Gas outlet |
| f Compartment for electrical components | |

All independant® Care Horizontal profiles are modular and allow a wide choice of configuration:



- With none, one or two lighting profiles
- With none, one or two DIN rail profiles

Lighting and rail profiles can be combined at the same time.

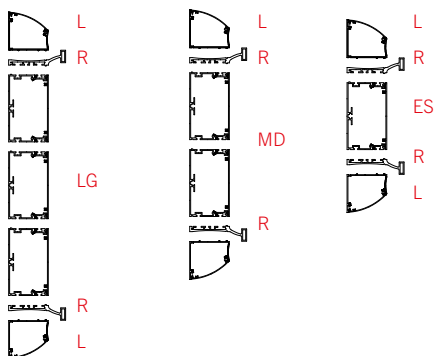


Vertical wall mounted


All vertical profiles are of a modular construction allowing a wide choice of configurations. Available upon request one or two stainless steel vertical poles onto which a number of dedicated accessories can be installed.

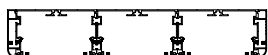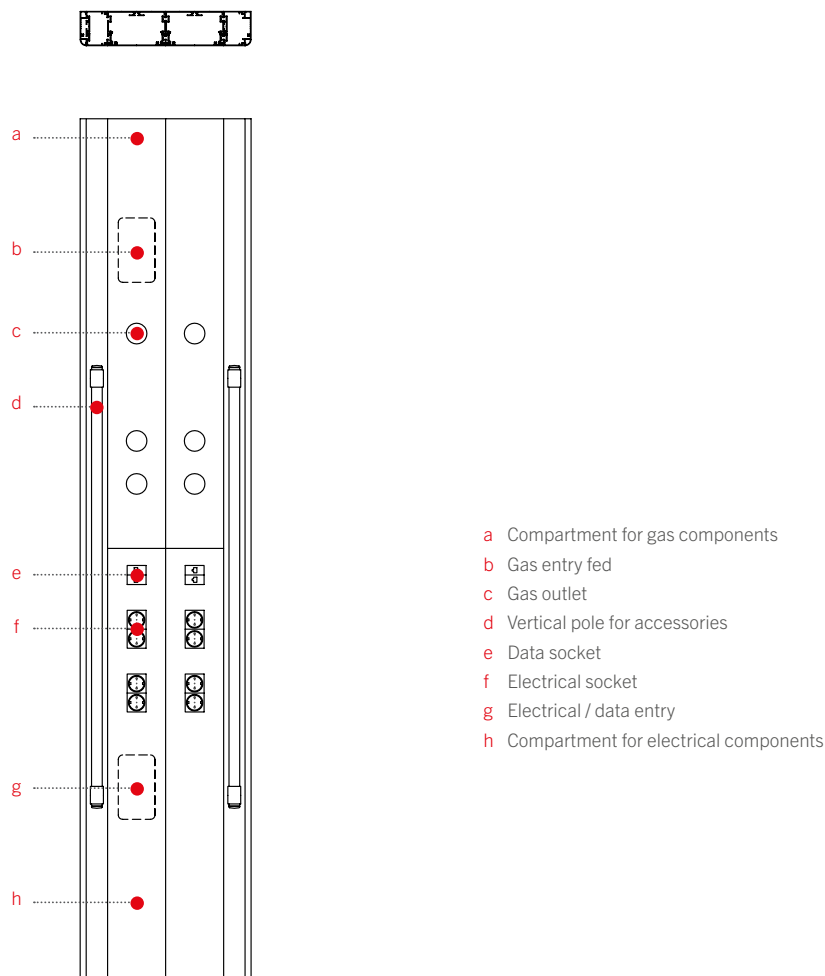

LG

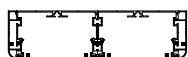


MD

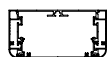

ES


This overview shows the various configurations available for each VERTICAL wall-mounted model. All VERTICAL profiles are of a modular construction allowing a wide choice of configurations:

- with one profile
- with two profiles
- with three profiles
- with multiple profiles (on demand)




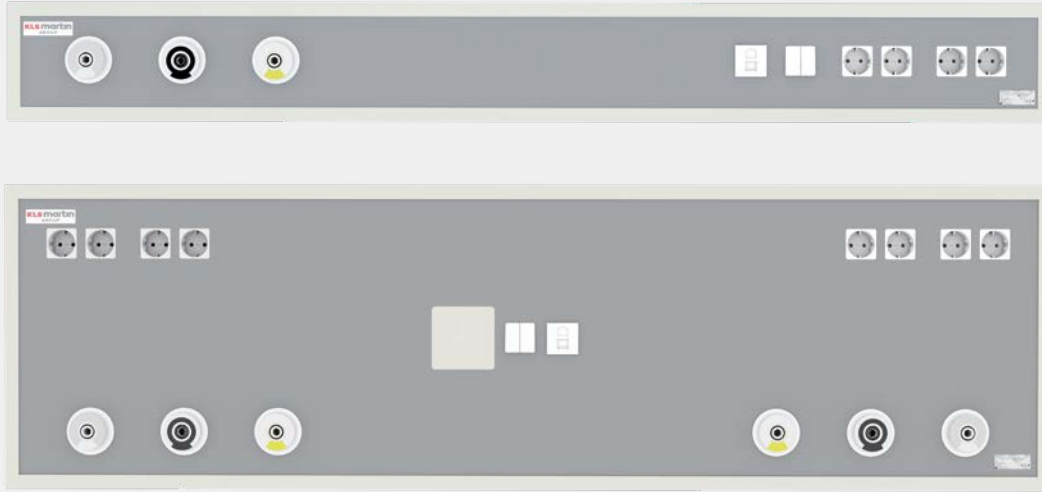


RECESSED MOUNTED

All models can be also recessed either in a concrete wall or any plasterboard. In addition to this, they can be completely integrated to excising HPL paneling. This is especially suitable where sleek but simple design needs to meet all medical requirements for standard specialty wards. The configuration can be either horizontal or vertical and the unit itself can provide single, double or multiple profile.

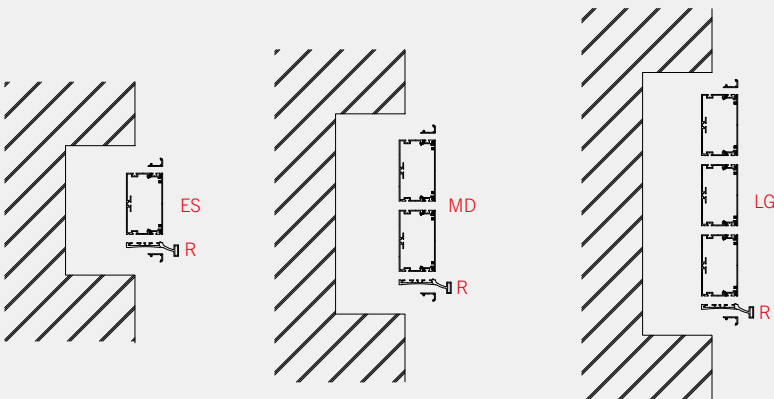


All configurations can be supplied with additional integrated DIN rail (horizontal mounting) or vertical pole (vertical mounting) for which a dedicated list of accessories is available.

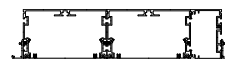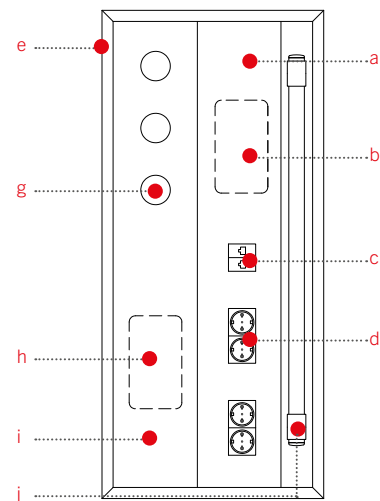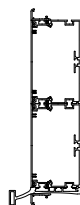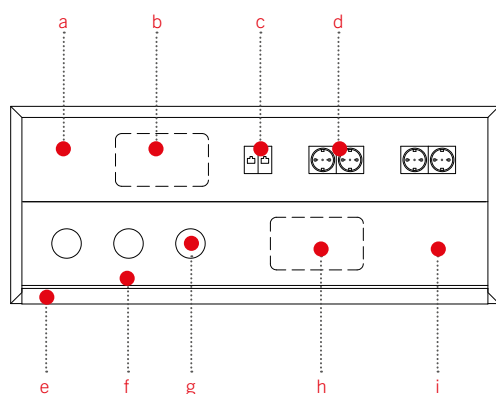


- a Compartment for electrical components
- b Electrical entry fed
- c Data socket
- d Electrical socket
- e Frame
- f DIN rail
- g Gas outlet
- h Gas entry fed
- i Compartment for gas components
- j Vertical pole for accessories

The recessed configuration gives a number of solutions that can fulfil any architectural and technical requirements. From single to two/multiple beds these units can be configured vertically or horizontally with one, two and even more profiles. Continuous twin solution is also available to serve two beds with one only solution.


Horizontal twin solution

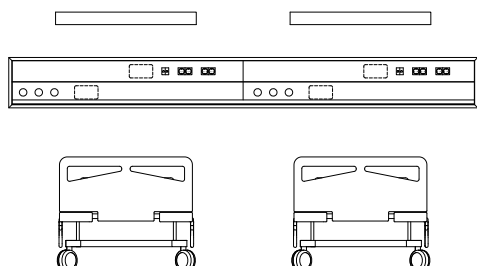

One-sided vertical solution

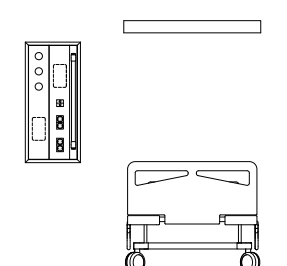

Double-sided vertical solution

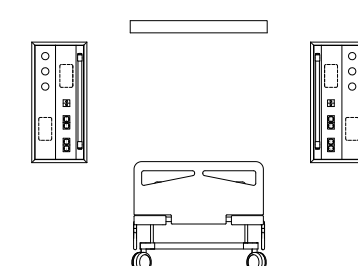

KLS Martin Group

KLS Martin Australia Pty Ltd.

Sydney · Australia
Tel. +61 2 9439 5316
australia@klsmartin.com

KLS Martin do Brasil Ltda.

São Paulo · Brazil
Tel. +55 11 3554 2299
brazil@klsmartin.com

KLS Martin Medical (Shanghai) International Trading Co., Ltd

Shanghai · China
Tel. +86 21 5820 6251
info@klsmartin.com

KLS Martin India Pvt Ltd.

Chennai · India
Tel. +91 44 66 442 300
india@klsmartin.com

KLS Martin Italia S.r.l.

Milan · Italy
Tel. +39 039 605 67 31
info@klsmartin.com

KLS Martin Japan K.K.

Tokyo · Japan
Tel. +81 3 3814 1431
info@klsmartin.com

KLS Martin SE Asia Sdn. Bhd.

Penang · Malaysia
Tel. +604 261 7060
malaysia@klsmartin.com

KLS Martin de México, S.A. de C.V.

Mexico City · Mexico
Tel. +52 55 7572 0944
mexico@klsmartin.com

KLS Martin Nederland B.V.

Huizen · Netherlands
Tel. +31 35 523 45 38
info@klsmartin.com

KLS Martin SE & Co. KG

Moscow · Russia
Tel. +7 499 792 76 19
russia@klsmartin.com

KLS Martin Taiwan Ltd.

Taipei · Taiwan
Tel. +886 2 2325 3169
taiwan@klsmartin.com

KLS Martin SE & Co. KG

Dubai · United Arab Emirates
Tel. +971 4 454 16 55
middleeast@klsmartin.com

KLS Martin UK Ltd.

Reading · United Kingdom
Tel. +44 118 467 1500
info.uk@klsmartin.com

KLS Martin LP

Jacksonville · Florida, USA
Tel. +1 904 641 77 46
usa@klsmartin.com

KLS Martin SE Asia Sdn. Bhd.

Hanoi · Vietnam
Tel. +49 7461 706-0
info@klsmartin.com

KLS Martin SE & Co. KG

A company of the KLS Martin Group

KLS Martin Platz 1 · 78532 Tuttlingen · Germany
PO Box 60 · 78501 Tuttlingen · Germany
Tel. +49 7461 706-0 · Fax +49 7461 706-193
info@klsmartin.com · www.klsmartin.com